


SKY pro guide – overview of drill application

Implant	Pilot drill	Twist drill	Final drill		Crestal drill	Implant mounter	
			Hard bone D1	Medium hard / soft bone D2, D3, D4		SKY	copaSKY
5.2 mm	4.0	PCDXLXD450	TDXL52D225	1. FD1L52D350 2. FD1L52D400	CRDXLXD400	n.a.	COPMLXD475
	4.5		TDXL52D225	1. FD1L52D350 2. FD1L52D400 3. FD1L52D450	CRDXLXD450	n.a.	COPMLXD475
8 mm	3.5 N	PCDXLXD450	1. TDXL52D225 2. TDXL08D225	1. FD1L52D350 2. FD1L08D350	CRDXLXD35n	SKYMLXD475	COPMLXD475
	4.0		1. TDXL52D225 2. TDXL08D225	1. FD1L52D350 2. FD1L52D400 3. FD1L08D400	CRDXLXD400	SKYMLXD475	COPMLXD475
	4.5		1. TDXL52D225 2. TDXL08D225	1. FD1L52D350 2. FD1L52D400 3. FD1L52D450 4. FD1L08D450	CRDXLXD450	SKYMLXD475	COPMLXD475
10 mm	3.5 N	PCDXLXD450	1. TDXL52D225 2. TDXL10D225	1. FD1L52D350 2. FD1L10D350	CRDXLXD35n	SKYMLXD475	COPMLXD475
	4.0		1. TDXL52D225 2. TDXL10D225	1. FD1L52D350 2. FD1L52D400 3. FD1L10D400	CRDXLXD400	SKYMLXD475	COPMLXD475
	4.5		1. TDXL52D225 2. TDXL10D225	1. FD1L52D350 2. FD1L52D400 3. FD1L52D450 4. FD1L10D450	CRDXLXD450	SKYMLXD475	COPMLXD475

Case study: * bSKY4514

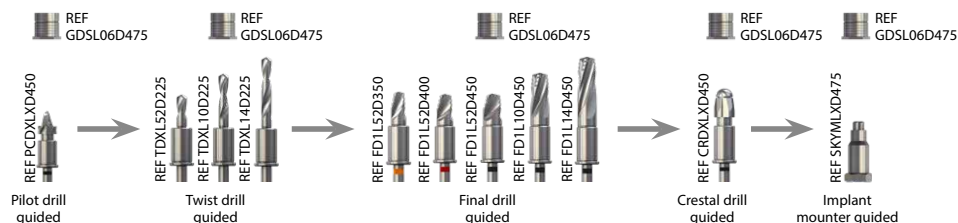

SKY pro guide – overview of drill application


Implant	Pilot drill	Twist drill	Final drill		Crestal drill	Implant mounter	
			Hard bone D1	Medium hard / soft bone D2, D3, D4		SKY	copaSKY
12 mm	3.5 N	PCDXLXD450	1. TDXL52D225 2. TDXL12D225	1. FD1L52D350 2. FD1L12D350	CRDXLXD35n	SKYMLXD475	COPMLXD475
	4.0		1. TDXL52D225 2. TDXL12D225	1. FD1L52D350 2. FD1L52D400 3. FD1L12D400	CRDXLXD400	SKYMLXD475	COPMLXD475
	4.5		1. TDXL52D225 2. TDXL12D225	1. FD1L52D350 2. FD1L52D400 3. FD1L52D450 4. FD1L12D450	CRDXLXD450	SKYMLXD475	COPMLXD475
14 mm	3.5 N	PCDXLXD450	1. TDXL52D225 2. TDXL10D225 3. TDXL14D225	1. FD1L52D350 2. FD1L10D350 3. FD1L14D350	CRDXLXD35n	SKYMLXD475	COPMLXD475
	4.0		1. TDXL52D225 2. TDXL10D225 3. TDXL14D225	1. FD1L52D350 2. FD1L52D400 3. FD1L10D400 4. FD1L14D400	CRDXLXD400	SKYMLXD475	COPMLXD475
	4.5 *		1. TDXL52D225 2. TDXL10D225 3. TDXL14D225	1. FD1L52D350 2. FD1L52D400 3. FD1L52D450 4. FD1L10D450 5. FD1L14D450	CRDXLXD450	SKYMLXD475	COPMLXD475
16 mm	3.5 N	PCDXLXD450	1. TDXL52D225 2. TDXL10D225 3. TDXL16D225	1. FD1L52D350 2. FD1L10D350 3. FD1L16D350	CRDXLXD35n	SKYMLXD475	n.a.
	4.0		1. TDXL52D225 2. TDXL10D225 3. TDXL16D225	1. FD1L52D350 2. FD1L52D400 3. FD1L10D400 4. FD1L16D400	CRDXLXD400	SKYMLXD475	n.a.




SKY
PRO GUIDE



Case study: copa4014

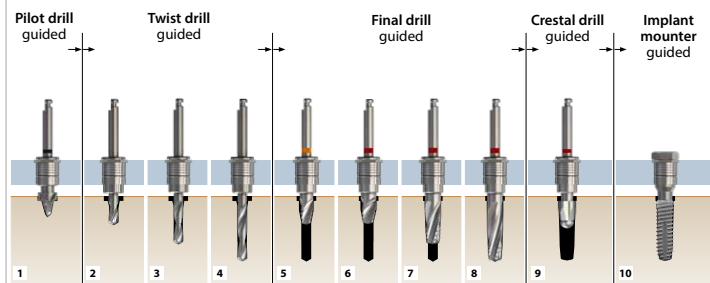

3.5 N
4.0
4.5

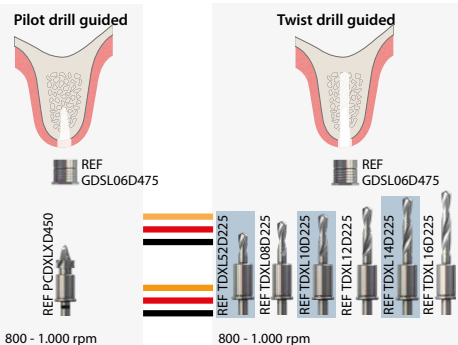

Implant length

8, 10, 12,
14, 16 mm

3.5 N

5, 2, 8, 10,
12, 14, 16 mm

4.0

5, 2, 8, 10,
12, 14, 16 mm

4.5

Implant mounter guided




Crestal drill guided





Hard bone
D1

Final drill guided

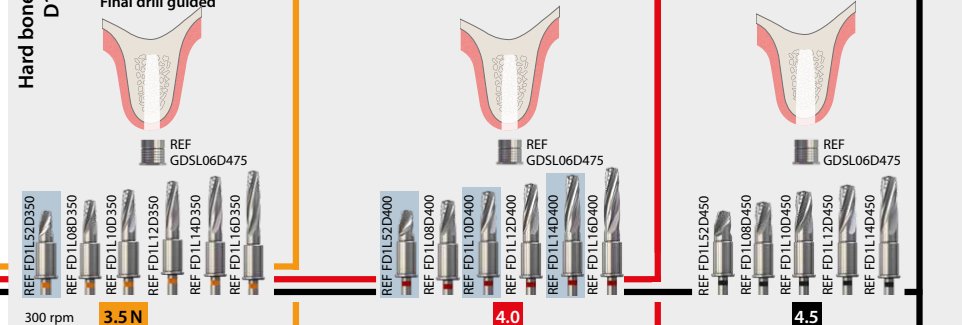

Medium hard /
soft bone
D2, D3, D4

Final drill guided

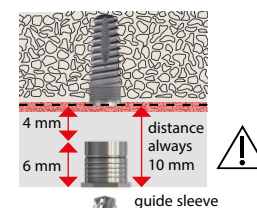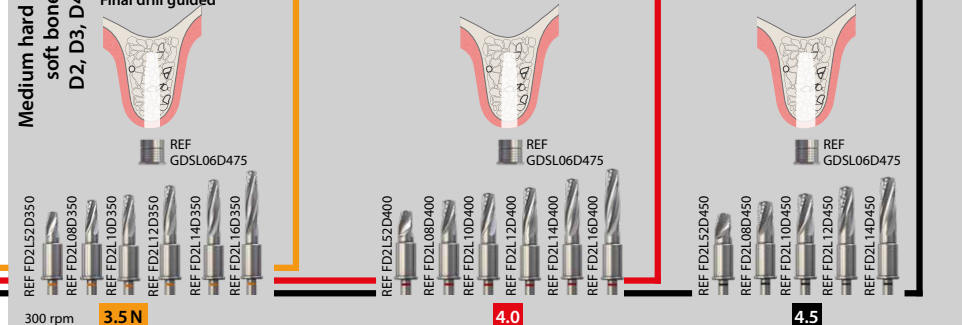

REF GDLSL06D475 REF GDLSL04D475

bredent
medical

bredent
group