


Milling unit BF 2

Precision guides ensure perfect milling and drilling results

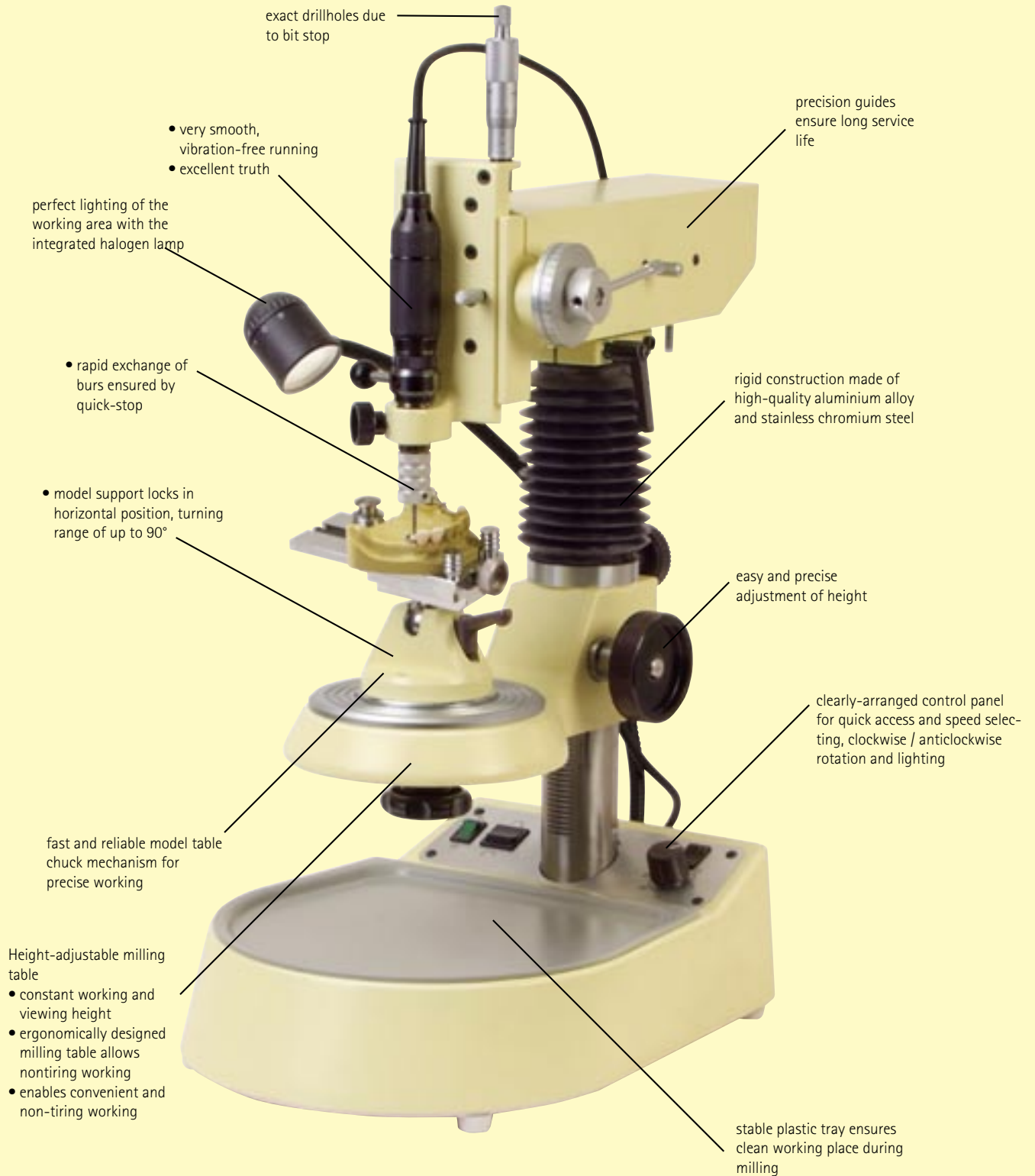


High-precision drilling and milling due to variably adjustable speed range from 0 to 30,000 rpm.

Milling unit BF 2


The milling unit that is simply perfect – even in its price/performance ratio



Excellent milling of wax due to variably adjustable speed.

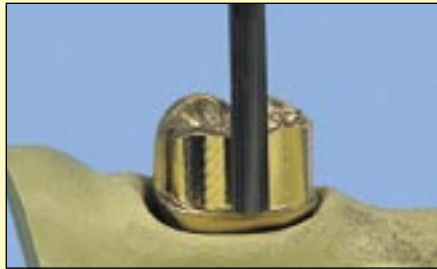

The bredent polishing bur allows to achieve shining metal surfaces.

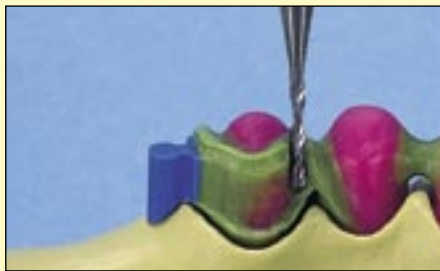

Interlocks and drillings in wax without clogging of the drillholes due to selection of low speed.



The high speed ensures safe and precise drilling, even in hard CoCr alloys.

Assortment

- 4 pieces
 - 1 Milling unit BF 2
 - 1 Handpiece BF 2
 - 1 Model support BF 1
 - 1 control unit BF 2
- REF 140 0098 0**



Marking the claps and surveying with a single unit – this is how time and money can be saved.



Placing the attachment, scraping wax surfaces parallel or conical – all this can be done with the BF 2.

Technical Data

Power supply	230 volts/ 50/60 Hz
Power rating	80 watts
Speed	0 – 30.000 rpm
Chuck	Ø 2.35 mm
Fuse	thermal overload protection
Torque	2.6 Ncm
Weight	18 kg
Width/Depth/Height	250 x 370 x 510 mm

Accessories

Chuck 2,35 mm	REF 730 0016 9
Chuck 3 mm	REF 730 0015 3
Tap handwheel	REF 330 0115 4
Model support BF 1	REF 730 0017 0
Milling base	REF 140 0089 3
Adapter airaqua turbine	
16 mm	REF 730 0018 4
18 mm (for BF1)	REF 730 0018 3
28,5 mm	REF 730 0018 5
Transfer device	
3 mm shaft	REF 360 0116 3
2,5 mm shaft	REF 360 0126 5



Model support BF 1



Milling base



Adapter airaqua turbine



Transfer device

Brenometer surveying system

