


PROCEDURA D' INSERIMENTO BIBLIOTECHE Attacchi BREDEnt in 3shape

1. Aprire Dental Manager e cliccare su “Control Panel”:

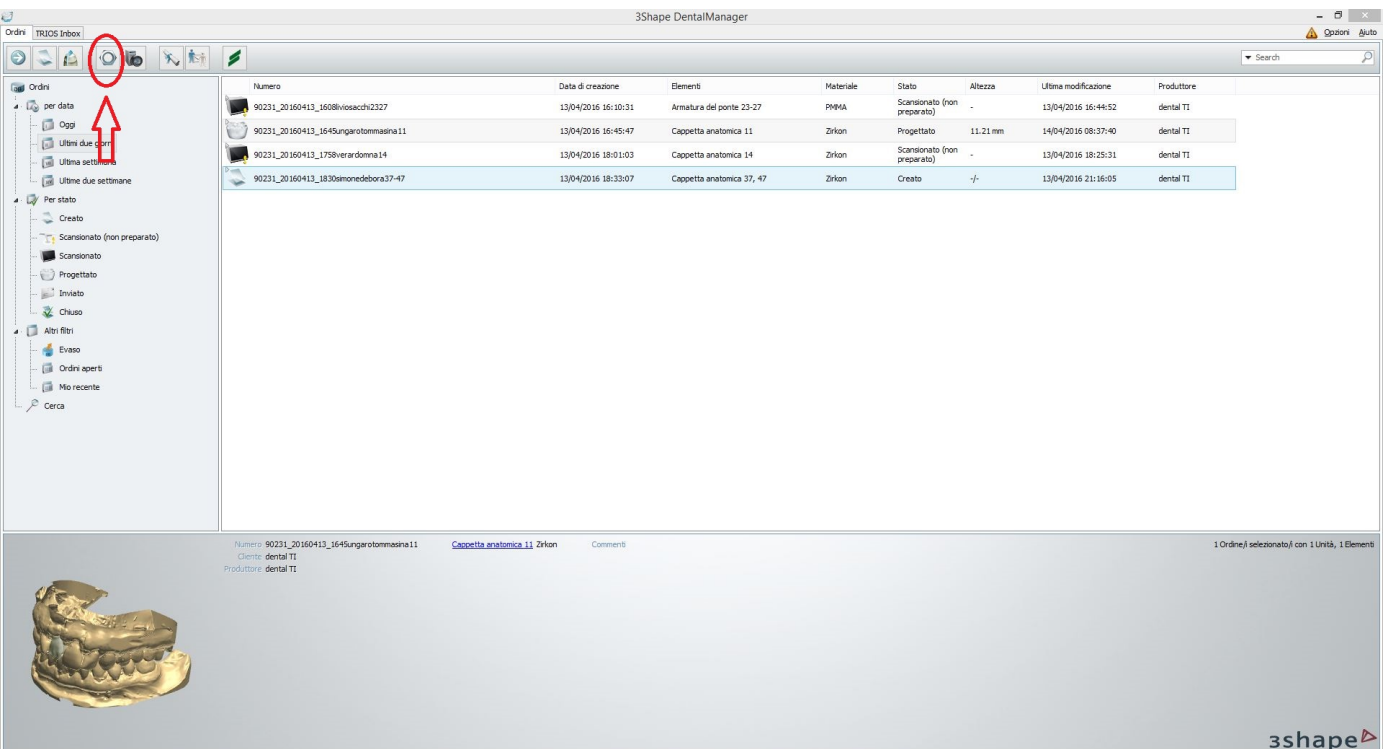


The screenshot shows the 3Shape DentalManager application window. The top toolbar contains several icons, with the 'Control Panel' icon (a gear) circled in red and an arrow pointing to it. The left sidebar shows a tree view of the application's structure. The main area displays a table of dental orders.

Numero	Data di creazione	Elementi	Materiale	Stato	Altezza	Ultima modificazione	Produttore
90231_20160413_1608lvioschi2327	13/04/2016 16:10:31	Armature del ponte 23-27	PMMA	Scansionato (non preparato)	-	13/04/2016 16:44:52	dental TI
90231_20160413_1645lungarotommasina11	13/04/2016 16:45:47	Cappetta anatomica 11	Zirkon	Progettato	11.21 mm	14/04/2016 08:37:40	dental TI
90231_20160413_1758verardonna14	13/04/2016 18:01:03	Cappetta anatomica 14	Zirkon	Scansionato (non preparato)	-	13/04/2016 18:25:31	dental TI
90231_20160413_1830simonedebora37-47	13/04/2016 18:33:07	Cappetta anatomica 37, 47	Zirkon	Creato	-/-	13/04/2016 21:16:05	dental TI

At the bottom left, there is a 3D model of a dental arch. The bottom right corner features the 3shape logo.

2. Cliccare su “Allegati”:

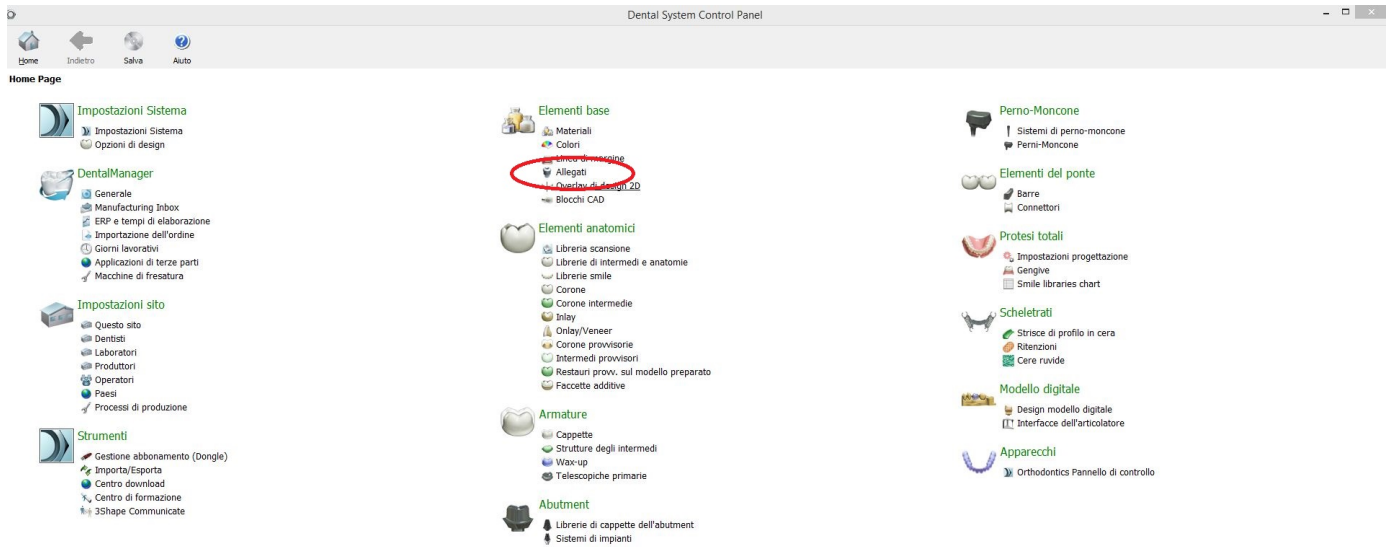

3. Cliccare su “Aggiungi” ed importare il file scaricato

